

## DATASHEET

**Product Code:** BRAPRPC7880 (PC) - BRAPRMO7880 (Monitor)

**Product Description:** ARTICULATED MONITOR AND PC ARMS

**Declaration of conformity:**

- CE
- Medical Device Class I : 1198796 - 1198819
- CND Code: V0899



### MISSION:

B&D articulated arms have a maximum capacity of 8,5 kg and are available in standard or with extension (380 mm) down the monitor. All the arms are equipped with VESA 7.5 and 10 mm with single vertical handle; this allows you to place it in a central location or in a position of + / - 25 mm. The connecting pivot has a diameter of 35mm or can be produced according to customer design and requests. Monitor arm version is equipped with cable 5 meter RCA - VGA for power and speaker, this allocation may be modified on request.




PC arm version is equipped with 5-meter network and power cable (diameter = 2.5): also this allocation may be amended at the request of the customer. In addition to USB cables, HDMI, S-VHS, other types of cables are available on customer request. All the movements are clutch, and this characteristics allows to adapt all the versions of our arms to any kind of devices of different weights. The same arm can be used with the pin facing upwards (ceiling) or down (to the floor), rotating towards its positioning on the central hub.

The fixed arm can be prepared on demand up to a length of 1200 mm.

**USE:**

- Monitor / Pc for patient Entertainment.
- Operating rooms, endoscopy Monitor or cameras systems.

**MOVEMENTS AND MEASURES**

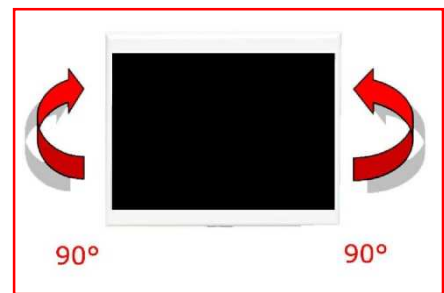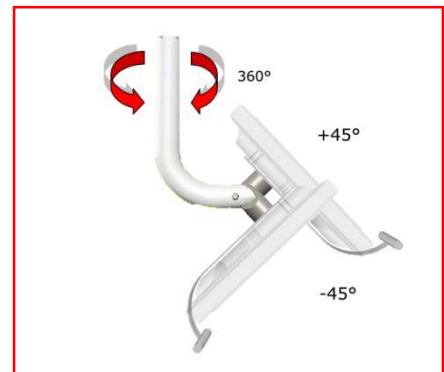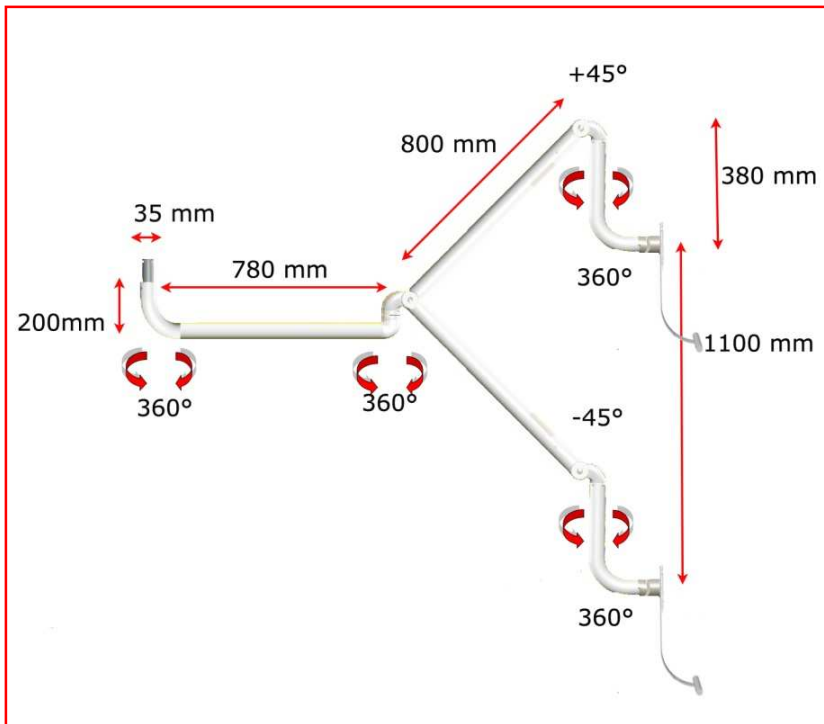

**FEATURES**

Arm Overall length	1780 mm
Weight bearable (Factory setting)	0,5 to 8,8 KG
Tilt Height (Vertical Movement)	1100 mm
Integrated cable	Monitor: Power - Lan PC : RCA – Power –VGA- Speaker
Handle with Vesa	75 -100 mm
Pivot:	Diameter: 35 mm
Mounting (Factory setting)	Ceiling or Wall / Floor
Total weight	7,5 KG
Standard color	RAL 9010
Customizations:	Cables - Color - Length fixed arm - length articulated arm - Limitation rotations -

Ver: 1.2