


PT MASTER 3's main specifications:

Frequency: 28000Hz	Tips hardness: 30HRC
Vibration amplitude: 20~60um	Treatment depth:14mm



Treatment effect:




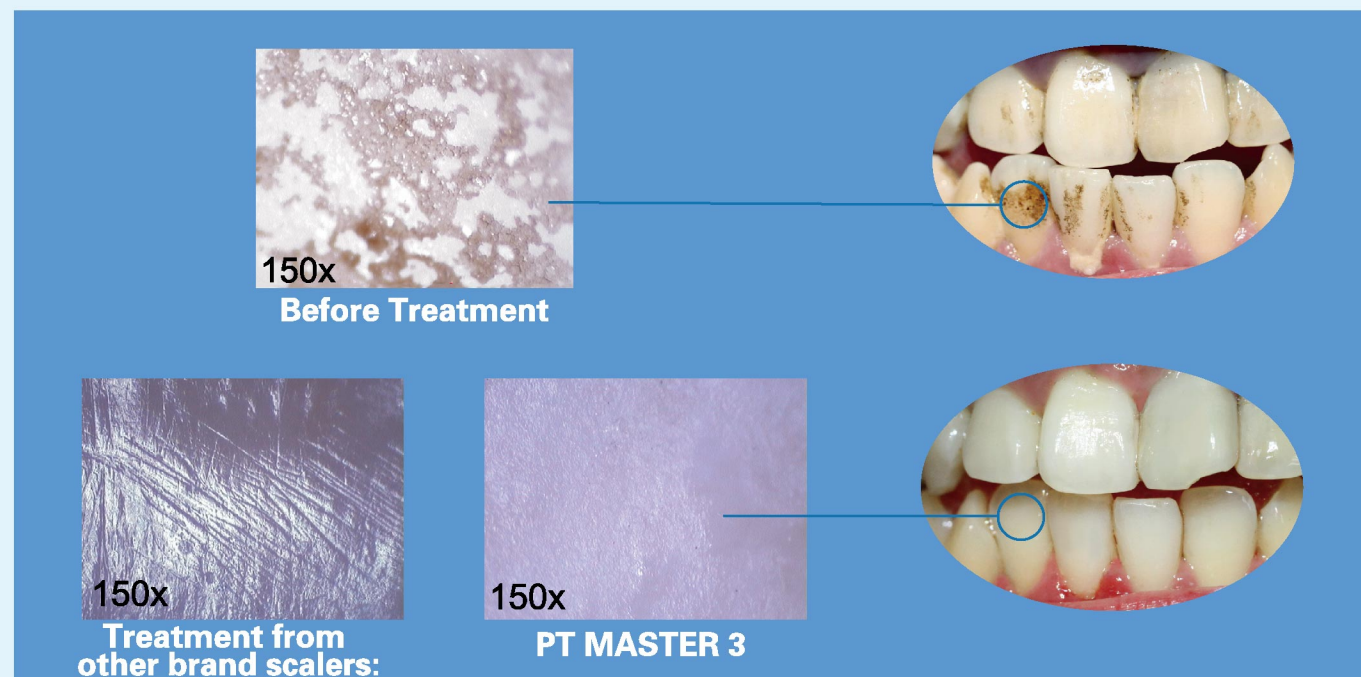




Less injury to teeth

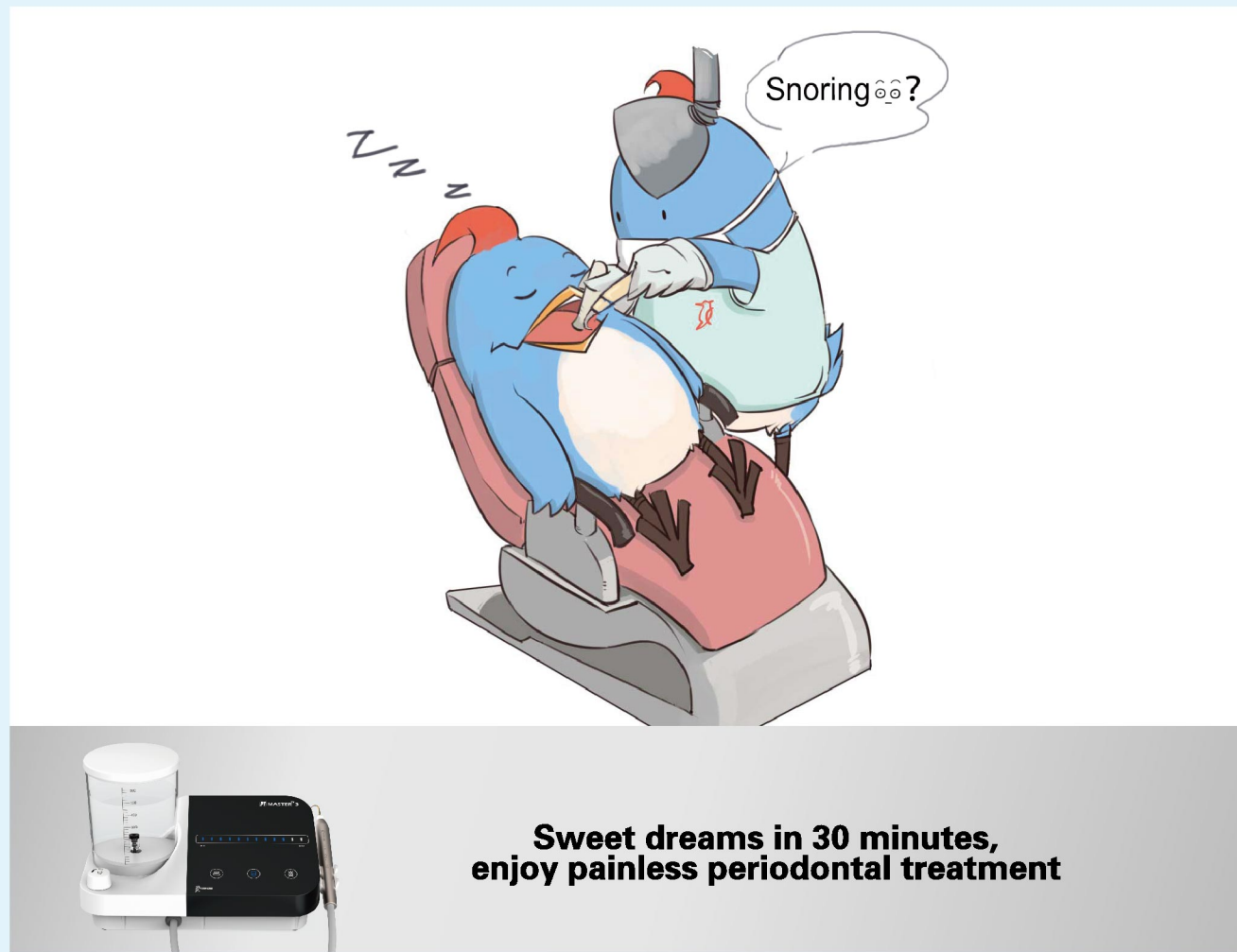

Contrast between PT MASTER 3 and Magnetostrictive device

	PT MASTER 3	Magnetostrictive device
Permanence of tips	Not easily broken 	Easily broken 
Price of tip	200-300Yuan/pc, (33-49 USD/pc)	2000-3000Yuan/pc, (330-490 USD/pc)
Safety and health of handpiece	Autoclavable and safe handpiece, avoid cross infection	Not fit to autoclave



-  **Periodontal treatment expert** | Circular vibration orbit, scaling and polishing together
-  **Titanium scaling tips** | No injury to cementum
-  **Precise location** | LED handpiece, better visibility, say no to accidental injury
-  **Low subsequent using cost** | PT MASTER3' s tough and preferential tip VS fragile highly cost magnetostrictive tip
-  **Painless treatment** | 20 μmvibrating amplitude, enjoy the painless treatment

Guilin Woodpecker Medical Instrument Co., Ltd.
 Add: Information Industrial Park, National High-Tech Zone, Guilin, Guangxi, 541004 P. R. China
 Tel: Europe Sales Dept.:+86-773-5873196, +86-773-2125222
 North America, South America & Oceania Sales Dept.: +86-773-5873198, +86-773-2125123
 Asia & Africa Sales Dept.:+86-773-5855350, +86-773-2125896
 Fax: +86-773-5822450
 E-mail: woodpecker@glwoodpecker.com, sales@glwoodpecker.com
 Web: http://www.glwoodpecker.com
 ZMN/WI-10-42 1.0 Edition

Distributed by:



Scan and Login website for more information



GLOBAL INITIAL CREATION

PIEZOELECTRIC CERAMICS PAINLESS PERIODONTAL TREATMENT DEVICE



1 Moveable water bottle

- Available to Hydrogen peroxide, sodium hypochlorite, chlorhexidine and special liquid
- Moveable water bottle, easier to change water
- Clear 20ml marks, convenient to match medical liquid

2 Water-control switch

3 Smart touching technique system

- Clear control panel
- Ultra-sensitive
- Waterproof design

4 Smart power indicator

5 Automatic water supply



Patented appearance, counterfeiting not allowed!
Patent No.: CN 201430379561.8

6 Periodontal treatment tips

- Special design for periodontal treatment, convenient for clinic operation
- Titanium scaling tips, no injury to cementum

7 Detachable LED handpiece

- Surprising efficiency
- Unique assemble technology, durable using life, far ahead of the other similar products

8 Cable

- High flexibility, durable using life

Circular vibration orbit of tip, special design for periodontal treatment




Perfect design for tips, 360° scaling



Usage of tip



P50R right-angled

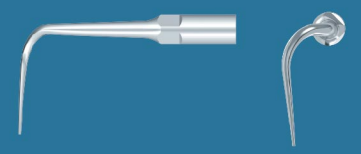

Used for subgingival scaling, periodontal pocket probing and root planning in the first and third quadrant.


P50L left-angled



Used for subgingival scaling, periodontal pocket probing and root planning in the second and fourth quadrant.



Easy to change water

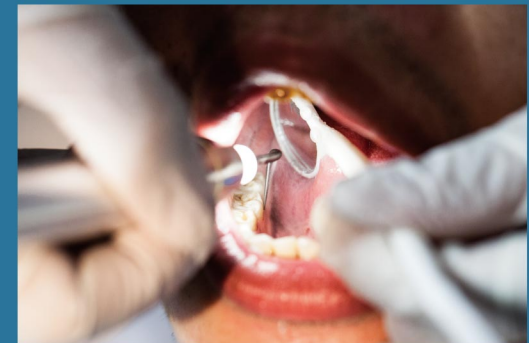

LED handpiece, better visibility



Ultra-sensitive control, waterproof design



Equipped with two professional disinfection suits